

Europäisches  
Patentamt

European Patent  
Office


Office européen  
des brevets

EP04/502029

REC'D 13 OCT 2004

WIPO

PCT

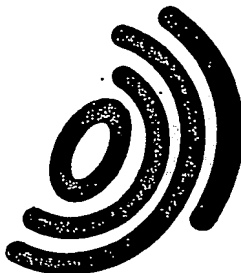

## Bescheinigung

## Certificate

## Attestation

Die angehefteten Unterlagen stimmen mit der ursprünglich eingereichten Fassung der auf dem nächsten Blatt bezeichneten internationalen Patentanmeldung überein.

The attached documents are exact copies of the international patent application described on the following page, as originally filed.

Les documents fixés à cette attestation sont conformes à la version initialement déposée de la demande de brevet international spécifiée à la page suivante.

Den Haag, den  
The Hague,  
La Haye, le

12. 10. 2004

Der Präsident des Europäischen Patentamts  
Im Auftrag  
For the President of the European Patent Office  
Le Président de l'Office européen des brevets  
p.o.

H.A.M.W. ter Haar

**PRIORITY  
DOCUMENT**  
SUBMITTED OR TRANSMITTED IN  
COMPLIANCE WITH RULE 17.1(a) OR (b)

Patentanmeldung Nr.  
Patent application no.  
Demande de brevet n°

PCT/EP 03/10263

**blatt 2 der Bescheinigung  
heet 2 of the certificate  
age 2 de l'attestation**

---

**Anmeldung Nr.:  
Application no.:  
Démarche n°:**

PCT/EP 03/10263

**Anmelder:  
Applicant(s):  
Démarcheur(s):**

1. JANSSEN PHARMACEUTICA N.V. - Beerse, Belgium  
2. SMANS, Karine, Alfonsine, Astrid - Beerse, Belgium (US only)

**Bezeichnung der Erfindung:  
Title of the invention:  
Titre de l'invention:**

CHIMERIC GABAB RECEPTOR

**Anmeldetag:  
Date of filing:  
Date de dépôt:**

12 September 2003 (12.09.2003)

**Anspruch genommene Priorität(en)  
Priority(ies) claimed  
Priorité(s) revendiquée(s)**

<b>Land: State: Pays:</b>	<b>Tag: Date: Date:</b>	<b>Aktenzeichen: File no. Numéro de dépôt:</b>
-----------------------------------	---------------------------------	--

**Benennung von Vertragsstaaten : Siehe Formblatt PCT/RO/101 (beigefügt)  
Designation of contracting states : See Form PCT/RO/101 (enclosed)  
Désignation d'états contractants : Voir Formulaire PCT/RO/101 (ci-joint)**

**Bemerkungen:  
Remarks:  
Remarques:**

## PCT REQUEST

Original (for SUBMISSION) - printed on 11.09.2003 10:22:38 AM

IV-1	<b>Agent or common representative; or address for correspondence</b> The person identified below is hereby/has been appointed to act on behalf of the applicant(s) before the competent International Authorities as:	<b>common representative</b>
IV-1-1	Name	<b>JANSSEN PHARMACEUTICA N.V.</b>
IV-1-2	Address:	<b>Patent Department Turnhoutseweg 30 B-2340 Beerse Belgium</b>
IV-1-3	Telephone No.	<b>+ 32 14 60 21 86</b>
IV-1-4	Facsimile No.	<b>+ 32 14 60 54 91</b>
IV-1-5	e-mail	<b>patents@janbe.jnj.com</b>
V	<b>Designation of States</b>	
V-1	<b>Regional Patent</b> (other kinds of protection or treatment, if any, are specified between parentheses after the designation(s) concerned)	<b>EP: AT BE BG CH&amp;LI CY CZ DE DK EE ES FI FR GB GR HU IE IT LU MC NL PT RO SE SI SK TR and any other State which is a Contracting State of the European Patent Convention and of the PCT</b>
V-2	<b>National Patent</b> (other kinds of protection or treatment, if any, are specified between parentheses after the designation(s) concerned)	<b>US</b>
V-5	<b>Precautionary Designation Statement</b>  In addition to the designations made under items V-1, V-2 and V-3, the applicant also makes under Rule 4.9(b) all designations which would be permitted under the PCT except any designation(s) of the State(s) indicated under item V-6 below. The applicant declares that those additional designations are subject to confirmation and that any designation which is not confirmed before the expiration of 15 months from the priority date is to be regarded as withdrawn by the applicant at the expiration of that time limit.	
V-6	<b>Exclusion(s) from precautionary designations</b>	<b>NONE</b>
VI	<b>Priority claim</b>	<b>NONE</b>
VII-1	<b>International Searching Authority Chosen</b>	<b>European Patent Office (EPO) (ISA/EP)</b>

## CHIMERIC GABA<sub>B</sub> RECEPTOR

The present invention provides a novel method to identify substances that are agonists of GABA<sub>B</sub> receptors, using a <sup>3</sup>H-GABA binding assay in recombinant GABA<sub>B</sub>R1a/R2 receptor expressing cells.

### BACKGROUND OF THE INVENTION

GABA (γ-amino-butyric acid) is the most widely distributed amino acid inhibitory neurotransmitter in the central nervous system (CNS) activating two distinct families of receptors; the ionotropic GABA<sub>A</sub> and GABA<sub>C</sub> receptors for fast synaptic transmissions, and the metabotropic GABA<sub>B</sub> receptors governing a slower synaptic transmission.

GABA<sub>B</sub> receptors are members of the superfamily of seven transmembrane G-protein coupled receptors that are coupled to neuronal K<sup>+</sup> or Ca<sup>2+</sup> channels. Presynaptic GABA<sub>B</sub> receptor activation has generally been reported to result in the inhibition of Ca<sup>2+</sup> conductance, leading to a decrease in the evoked release of neurotransmitters. Post-synaptically the major effect of GABA<sub>B</sub> receptor activation is to open potassium channels, to generate post-synaptic inhibitory potentials.

The expression of GABA<sub>B</sub> receptors is widely distributed in the mammalian neuronal axis, with particularly high levels in the molecular layer of the cerebellum, interpeduncular nucleus, frontal cortex, olfactory nuclei, thalamic nuclei, temporal cortex, raphe magnus and spinal cord. GABA<sub>B</sub> receptors are also present in the peripheral nervous system, both on sensory nerves and on parasympathetic nerves. Their ability to modulate these nerves give them potential as targets in disorders of the lung, GI tract and bladder (Belley et al., 1999, *Biorg. Med. Chem.* 7:2697-2704).

A large number of pharmacological activities have been attributed to GABA<sub>B</sub> receptor activation, such as for example, analgesia, hypothermia, catatonia, hypotension, reduction of memory consolidation and retention, and stimulation of insulin, growth hormone and glucagon release (see Bowery, 1989, *Trends Pharmacol. Sci.* 10:401-407 for a review). It is well accepted that GABA<sub>B</sub> receptor agonists and antagonists are pharmacologically useful in indications such as stiff man syndrome, gastroesophageal reflux, neuropathic pain, incontinence and treatment of cough and cocaine addiction. For example, the GABA<sub>B</sub> receptor agonist baclofen has been shown to reduce transient lower esophageal sphincter relaxations (TLESR) and is accordingly useful in the treatment of reflux as most episodes of reflux occur during TLESR. However, the current GABA<sub>B</sub> receptor agonists, such as baclofen, are relatively non-selective and show a variety of undesirable behavioural actions such as sedation and

respiratory depression. It would be desirable to develop more GABA<sub>B</sub> receptor agonists with an improved selectivity and less of the aforementioned undesirable effects.

5 Current methods of drug discovery generally involve assessing the biological activity of tens or hundreds of thousands of compounds in order to identify a small number of those compounds having a desired activity against a particular target, i.e. High Throughput Screening (HTS). In a typical HTS related screen format, assays are performed in multi-well microplates, such as 96, 384 or 1536 well plates, putting certain constraints to the setup of the assay to be performed including the availability of  
10 the source materials (i.e membrane preparations of cells expressing the recombinant GABA<sub>B</sub> receptor). HTS related screens are preferably performed at room temperature with a single measurement for each of the compounds tested in the assay, requiring short cycle times, with a reproducible and reliable output.

15 Present *in vitro* screens to identify compounds as agonists of the GABA<sub>B</sub> receptor, either rely on natural, less abundant resources such as binding assays in rat brain membranes or consist of functional screening assays, such as for example Ca<sup>2+</sup> responses, c-AMP responses and effects on Ca<sup>2+</sup> and K<sup>+</sup> channels performed in cells expressing a recombinant GABA<sub>B</sub> receptor. In some of these functional assays the GABA<sub>B</sub> receptors may be co-expressed with G-proteins, e.g. Gα16 or Gq15 or the  
20 chimeric G-protein G αq-z5, increasing G-protein coupling (Bräuner-Osborne & Krogsgaard-Larsen; 1999, Br. J. Pharmacol. 128:1370-1374). However, a GABA<sub>B</sub> agonist binding assay that would further reduce the HTS cycle time and the resources for biochemicals such as recombinant proteins, is currently unavailable.

25 The present invention describes the development of a Chinese Hamster Ovary (CHO) cell line co-expressing the human GABA<sub>B</sub> receptor subunits GABA<sub>B</sub>R1a and GABA<sub>B</sub>R2, which were surprisingly found to demonstrate agonist binding in radioligand binding experiments. In addition, the present inventors demonstrated that the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line has one high affinity and one low affinity agonist binding site in the recombinant expressed GABA<sub>B</sub> receptor. Hence the  
30 hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line provided by the present invention not only allows compound screening, but also provides a useful tool to characterize the nature of the compound –receptor interaction.

## SUMMARY OF THE INVENTION

The present invention provides an isolated GABA<sub>B</sub> receptor protein comprising at least one GABA<sub>B</sub>R1a subunit and at least one GABA<sub>B</sub>R2 subunit, characterized in that said GABA<sub>B</sub> receptor has one high affinity agonist binding site and one low affinity agonist binding site. In particular the isolated recombinant GABA<sub>B</sub> receptor protein expressed by the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone on August 22, 2003 with the accession number **LMBP XXXX**. It is thus an object of the present invention to provide the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone on August 22, 2003 with the accession number **LMBP XXXX**.

The invention also provides the use of the aforementioned cell line in a method to identify GABA<sub>B</sub> receptor agonists using a functional or a binding assay. In particular in a radioligand-binding assay comprising the use of radiolabeled agonists such as for example <sup>3</sup>H-GABA or <sup>3</sup>H-baclofen.

The invention further provides a method to identify GABA<sub>B</sub> receptor agonists, comprising contacting the aforementioned cell line with a test compound and measuring the binding of said test compound to the GABA<sub>B</sub> receptor. In particular the method consists of a radioligand binding assay, comprising exposing the aforementioned cells to a labelled agonist of GABA<sub>B</sub> in the presence and absence of the test compound and measure the binding of the labelled ligand to the cells according to the invention, where if the amount of binding of the labelled ligand is less in the presence of the test compound, then the compound is a potential agonist of the GABA<sub>B</sub> receptor.

It is also an object of the present invention to provide a method to identify a high affinity GABA<sub>B</sub> receptor agonist, said method comprising contacting the aforementioned cells with the radiolabeled agonist selected from the group consisting of GABA, baclofen and 3-aminopropylphosphinic acid (3-APPA a.k.a APMPA), in the presence and absence of the test compound and measure the binding of the labelled ligand to the cells according to the invention, where if the amount of binding of the labelled ligand to the high affinity binding site is less in the presence of the test compound, then the compound is a potential high affinity agonist of the GABA<sub>B</sub> receptor.

Alternatively, the aforementioned binding assays are performed on cellular extracts, in particular cellular membrane preparations of the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated

Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone on August 22, 2003 with the accession number **LMBP XXXX**.

5 In another embodiment the present invention provides a method to identify a GABA<sub>B</sub> receptor agonist, said method comprising contacting the aforementioned cell line with a compound to be tested and determine whether the compound activates a GABA<sub>B</sub> receptor functional response in said cells. In particular the functional response consists of modulation of the activity of ion channels or of intracellular messengers as explained hereinafter.

10 This and further aspects of the present invention will be discussed in more detail hereinafter.

#### BRIEF DESCRIPTION OF THE DRAWING

15 Figure 1 GTP $\gamma$ 35S-binding upon stimulation of membranes by GABA expressed as the percentage of maximal GABA stimulation, in the presence and absence of the positive allosteric modulator CGP7930>

20 Figure 2 Displacement of <sup>3</sup>H-GABA by agonists (baclofen, GABA & APMPA) and antagonists (SCH50911 & CGP54626)

Figure 3 Reproducible agonist IC<sub>50</sub> values (n=5) independent of membrane preparations.

25

#### DETAILED DESCRIPTION

For the purposes of describing the present invention: GABA<sub>B</sub>R1a or h GABA<sub>B</sub>R1a as used herein refers to the human GABA<sub>B</sub> receptor subunit known as GABA<sub>B</sub>R1a in Kaupmann et al, 1998, Proc. Natl. Acad. Sci. USA 95:14991-14996, the amino acid sequence (SEQ ID No.:2) of which can be found at GenBank Accession no. AJ225028, as well as to its mammalian orthologs. GABA<sub>B</sub>R1a also refers to other GABA<sub>B</sub> receptor subunits that have minor changes in amino acid sequence from those described hereinbefore, provided those other GABA<sub>B</sub> receptor subunits have substantially the same biological activity as the subunits described hereinbefore. A GABA<sub>B</sub>R1a subunit

30

35

has substantially the same biological activity if it has an amino acid sequence that is at least 80% identical to, preferably at least 95% identical to, more preferably at least 97% identical to, and most preferably at least 99% identical to SEQ ID No.: 2 and has a Kd or EC50 for GABA, GABA<sub>B</sub> receptor agonists such as for example baclofen and gabapentin or GABA<sub>B</sub> receptor antagonists such as for example CGP54626A, SCH 50911, saclofen and phaclofen, that is no more than 5-fold greater than the Kd or EC50 of a native GABA<sub>B</sub> receptor for GABA or the same GABA<sub>B</sub> receptor agonist or GABA<sub>B</sub> receptor antagonist.

10 GABA<sub>B</sub>R2 as used herein refers to the human GABA<sub>B</sub> receptor subunit known as GABA<sub>B</sub>R2 in White et al., 1998, Nature 396:679-682, the amino acid sequence (Seq Id NO.: 4) of which can be found at GenBank accession no. AF058795 as well as to its mammalian orthologs. GABA<sub>B</sub>R2 also refers to other GABA<sub>B</sub> receptor subunits that have minor changes in amino acid sequence from those described hereinbefore,  
15 provided those other GABA<sub>B</sub> receptor subunits have substantially the same biological activity as the subunits described hereinbefore. A GABA<sub>B</sub>R2 subunit has substantially the same biological activity if it has an amino acid sequence that is at least 80% identical to, preferably at least 95% identical to, more preferably at least 97% identical to, and most preferably at least 99% identical to SEQ ID No.: 4 and has in combination  
20 with a GABA<sub>B</sub>R1 subunit a Kd or EC50 for GABA, GABA<sub>B</sub> receptor agonists such as for example baclofen and gabapentin or GABA<sub>B</sub> receptor antagonists such as for example CGP54626A, SCH 50911, saclofen and phaclofen, that is no more than 5-fold greater than the Kd or EC50 of a native GABA<sub>B</sub> receptor for GABA or the same GABA<sub>B</sub> receptor agonist or GABA<sub>B</sub> receptor antagonist.

25 The Kd and EC50 values of the native GABA<sub>B</sub> receptor is determined using the methods known to a person skilled in the art, in particular using competition binding studies on tissue preparations such as for example described in Cross & Horton, 1987 Eur.J.Pharmacol. 141(1): 159-162. Briefly, crude synaptic membranes are prepared by  
30 homogenisation of whole brain, centrifugation (30 000 xg, 20 min.) and extensive washing. Total binding is measured by incubation of the membranes with <sup>3</sup>H-GABA or <sup>3</sup>H-baclofen, while non-specific binding is measured in the presence of 100 μM baclofen. Upon removal of unbound ligand by filtration, filters are counted in a β-counter or a Topcount Harvester (Packard). For competition experiments the binding  
35 occurs in the presence of increasing concentration of unlabeled compound.



It is thus an object of the present invention to provide an isolated GABA<sub>B</sub> receptor protein formed by at least one GABA<sub>B</sub>R1a and at least one GABA<sub>B</sub>R2 subunit further characterized in that said isolated GABA<sub>B</sub> has both a high and a low affinity agonist binding site. In a further embodiment this isolated GABA<sub>B</sub> receptor is a functional GABA<sub>B</sub> receptor expressed by a cell, wherein said cell does not normally express the GABA<sub>B</sub> receptor. Suitable cells which are commercially available, include but are not limited to L-cells, HEK-293 cells, COS cells, CHO cells, HeLa cells and MRC cells, in particular CHO cells wherein the GABA<sub>B</sub> receptor protein comprises at least one GABA<sub>B</sub>R1a subunit encoded by the oligonucleotide sequence consisting of SEQ ID No.1 and at least one GABA<sub>B</sub>R2 subunit encoded by the oligonucleotide sequence consisting of SEQ ID No.3. In a more particular embodiment the isolated GABA<sub>B</sub> receptor according to the invention, consists of the receptor protein expressed by the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone 20 on August 22, 2003 with accession number LMBP XXXX .

"Functional GABA<sub>B</sub> receptor" refers to a GABA<sub>B</sub> receptor formed by co-expression of GABA<sub>B</sub>R2 and GABA<sub>B</sub>R1a in a cell, wherein said cell does not normally express the GABA<sub>B</sub> receptor, most preferably resulting in a heterodimer of GABA<sub>B</sub>R2 and GABA<sub>B</sub>R1a, where the functional GABA<sub>B</sub> receptor mediates at least one functional response when exposed to the GABA<sub>B</sub> receptor agonist GABA. Examples of functional responses are : pigment aggregation in *Xenopus* melanophores, negative modulation of cAMP levels, coupling to inwardly rectifying potassium channels, mediation of late inhibitory postsynaptic potentials in neurons, increases in potassium conductance, decreases in calcium conductance, MAPKinase activation, extracellular pH acidification, and other functional responses typical of G-protein coupled receptors. One skilled in the art would be familiar with a variety of methods of measuring the functional responses of G-protein coupled receptors such as the GABA<sub>B</sub> receptor (see, e. g., Lerner, 1994, *Trends Neurosci.* 17: 142-146 [changes in pigment distribution in melanophore cells]; Yokomizo et al., 1997, *Nature* 387: 620-624 [changes in cAMP or calcium concentration; chemotaxis]; Howard et al., 1996, *Science* 273: 974-977 [changes in membrane currents in *Xenopus* oocytes]; McKee et al., 1997, *Mol. Endocrinol.* 11: 415-423 [changes in calcium concentration measured using the aequorin assay]; Offermanns & Simon, 1995, *J. Biol. Chem.* 270: 15175-15180 [changes in inositol phosphate levels]). Depending upon the cells in which heteromers of GABA<sub>B</sub>R1a and GABA<sub>B</sub>R2 are expressed, and thus the G-proteins with which the functional GABA<sub>B</sub> receptor thus formed is coupled, certain of such methods may be

appropriate for measuring the functional responses of such functional GABA<sub>B</sub> receptors. It is well within the competence of one skilled in the art to select the appropriate method of measuring functional responses for a given experimental system.

5           The term "compound", "test compound", "agent" or "candidate agent" as used herein can be any type of molecule, including for example, a peptide, a polynucleotide, or a small molecule that one wishes to examine for their activity as GABA<sub>B</sub> receptor agonist, and wherein said agent may provide a therapeutic advantage to the subject receiving it. The candidate agents can be administered to an individual by various  
10 routes, including, for example, orally or parenterally, such as intravenously, intramuscularly, subcutaneously, intraorbitally, intracapsularly, intraperitoneally, intrarectally, intracisternally or by passive or facilitated absorption through the skin, using for example a skin patch or transdermal iontophoresis, respectively. Furthermore the compound can be administered by injection, intubation or topically, the latter of  
15 which can be passive, for example, by direct application of an ointment, or active, for example, using a nasal spray or inhalant, in which case one component of the composition is an appropriate propellant. The route of administration of the compound will depend, in part, on the chemical structure of the compound. Peptides and polynucleotides, for example, are not particular useful when administered orally  
20 because they can be degraded in the digestive tract. However, methods for chemically modifying peptides, for example rendering them less susceptible to degradation are well known and include for example, the use of D-amino acids, the use of domains based on peptidomimetics, or the use of a peptoid such as a vinylogous peptoid.

          The agent used in the screening method may be used in a pharmaceutically acceptable carrier. See, e.g., *Remington's Pharmaceutical Sciences*, latest edition, by  
25 E.W. Martin Mack Pub. Co., Easton, PA, which discloses typical carriers and conventional methods of preparing pharmaceutical compositions that may be used in conjunction with the preparation of formulations of the agents and which is incorporated by reference herein.

30

### *Cells*

          As already outlined above, the present invention provides a cell line stably transfected with expression vectors that direct the expression of the GABA<sub>B</sub> receptor subunits GABA<sub>B</sub>R1a and GABA<sub>B</sub>R2 as defined hereinbefore. In particular CHO cells  
35 transfected with said expression vectors. Such expression vectors are routinely constructed in the art of molecular biology and may involve the use of plasmid DNA

and appropriate initiators, promoters, enhancers and other elements, which may be necessary, and which are positioned in the correct orientation, in order to allow for protein expression. Generally, any system or vector suitable to maintain, propagate or express polynucleotides to produce a polypeptide in a host may be used. The appropriate nucleotide sequence, i.e. the polynucleotide sequences encoding either the human GABA<sub>B</sub>R1a or GABA<sub>B</sub>R2 subunit as defined hereinbefore, may be inserted into an expression system by any of a variety of well-known and routine techniques such as for example those set forth in Current Protocols in Molecular Biology, Ausbel et al. eds., John Wiley & Sons, 1997.

In a particular embodiment the CHO cells according to the invention are cotransfected with the commercially available expression vectors pcDNA3.1 comprising the polynucleotide sequences encoding for human GABA<sub>B</sub>R1a (SEQ ID No.:1) and human GABA<sub>B</sub>R2 (SEQ ID No.: 3) respectively. More preferably the present invention provides a hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone 20 on August 22, 2003 with the accession number LMBP XXXX. This cell line is characterized in that the functional GABA<sub>B</sub> receptor in this CHO cell line has both a low and a high affinity binding site for GABA<sub>B</sub> receptor agonist. Using the cell line according to the invention, will not only allow compound screening, but also provides a useful tool for the characterization of the nature of the compound-receptor interaction, i.e. does it interact with the low or high affinity agonist binding site of the GABA<sub>B</sub> receptor.

For further details in relation to the preparation of nucleic acid constructs, mutagenesis, sequencing, introduction of DNA into cells and gene expression, and analysis of proteins, see for example, Molecular Cloning: a Laboratory Manual: 2nd edition, Sambrook et al., 1989, Cold Spring Harbor Laboratory Press.

### *Assays*

The present invention also provides an assay for a compound capable of interacting with the functional GABA<sub>B</sub> receptor of the present invention, which assay comprises: providing the GABA<sub>B</sub> receptor expressed by the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line of the present invention, contacting said receptor with a putative binding compound; and determining whether said compound is able to interact with said receptor.

In one embodiment of the assay, the receptor or subunits of the receptor may be employed in a binding assay. Binding assays may be competitive or non-competitive. Such an assay can accommodate the rapid screening of a large number of compounds to determine which compounds, if any, are capable of binding to the polypeptides.

Within this context, the present invention provides a method to identify whether a test compound binds to an isolated GABA<sub>B</sub> receptor protein of the present invention, and is thus a potential agonist or antagonist of the GABA<sub>B</sub> receptor, said method comprising;

a) contacting cells expressing a functional GABA<sub>B</sub> receptor, wherein such cells do not normally express the GABA<sub>B</sub> receptor, with the test compound in the presence and absence of a compound known to bind the GABA<sub>B</sub> receptor, and

b) determine the binding of the test compound to the GABA<sub>B</sub> receptor using the compound known to bind to the GABA<sub>B</sub> receptor as a reference.

Binding of the test compound or of the compound known to bind to the GABA<sub>B</sub> receptor, hereinafter also referred to as reference compound, is assessed using art-known methods for the study of protein-ligand interactions. For example, such binding can be measured by employing a labeled substance or reference compound. The test compound or reference compound can be labeled in any convenient manner known in the art, *e.g.* radioactively, fluorescently or enzymatically. In a particular embodiment of the aforementioned method, the compound known to bind to the GABA<sub>B</sub> receptor, a.k.a. the reference compound is detectably labeled, and said label is used to determine the binding of the test compound to the GABA<sub>B</sub> receptor. Said reference compound being labeled using a radiolabel, a fluorescent label or an enzymatic label, more preferably a radiolabel. In a more particular embodiment, the present invention provides a method to identify whether a test compound binds to an isolated GABA<sub>B</sub> receptor protein, said method comprising the use of a compound known to bind to the GABA<sub>B</sub> receptor, wherein said reference compound is selected from the group consisting of <sup>3</sup>H-GABA, <sup>3</sup>H-baclofen, <sup>3</sup>H-3-APPA, <sup>3</sup>H-CGP542626 and <sup>3</sup>H-SCH50911.

Subsequently, more detailed assays can be carried out with those compounds found to bind, to further determine whether such compounds act as agonists or antagonists of the polypeptides of the invention.

Thus, in a further embodiment the present invention provides a method to identify GABA<sub>B</sub> receptor agonists said method comprising,

- a) exposing cells expressing a functional GABA<sub>B</sub> receptor, wherein such cells do not normally express the GABA<sub>B</sub> receptor, to a labeled agonists of GABA<sub>B</sub> in the presence and absence of the test compound, and
- b) determine the binding of the labeled agonist to said cells,
- 5 where if the amount of binding of the labeled agonist is less in the presence of the test compound, then the compound is a potential agonist of the GABA<sub>B</sub> receptor. As already specified for the general binding assay above, the binding of the GABA<sub>B</sub> receptor agonists is assessed using art-known methods for the study of protein-ligand interactions. The label is generally selected from a radioactive label, a fluorescent label
- 10 or an enzymatic label, in particular a radiolabel wherein the agonist is selected from the group consisting of <sup>3</sup>H-GABA, <sup>3</sup>H-baclofen and <sup>3</sup>H-3-APPA.

Similarly, the present invention provides a method to identify GABA<sub>B</sub> receptor antagonists said method comprising,

- 15 a) exposing cells expressing a functional GABA<sub>B</sub> receptor, wherein said cells do not normally express the GABA<sub>B</sub> receptor, to a labeled antagonist of GABA<sub>B</sub> in the presence and absence of the test compound, and
- b) determine the binding of the labeled antagonist to said cells,
- 20 where if the amount of binding of the labeled antagonist is less in the presence of the test compound, then the compound is a potential antagonist of the GABA<sub>B</sub> receptor. As already specified for the general binding assay above, the binding of the GABA<sub>B</sub> receptor antagonists is assessed using art-known methods for the study of protein-ligand interactions. The label is generally selected from a radioactive label, a fluorescent label or an enzymatic label, in particular a radiolabel wherein the antagonist is selected from
- 25 the group consisting of <sup>3</sup>H-CGP542626 and <sup>3</sup>H-SCH50911.

In an alternative embodiment of the present invention, the aforementioned binding assays are performed on a cellular composition, i.e a cellular extract, a cell fraction or cell organelles comprising a GABA<sub>B</sub> receptor as defined hereinbefore. More in

30 particular, the aforementioned binding assays are performed on a cellular composition, i.e. a cellular extract, a cell fraction or cell organelles comprising a GABA<sub>B</sub> receptor as defined hereinbefore, wherein said cellular composition, i.e. cellular extract, cell fraction or cell organelles, is obtained from cells expressing a functional GABA<sub>B</sub> receptor, wherein said cells do not normally express the GABA<sub>B</sub> receptor. More

35 preferably, the cellular composition, i.e. cellular extract, cell fraction or cell organelles, is obtained from the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian

Coordinated Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone 20 on August 22, 2003 with the accession number **LMBP XXXX**.

It is accordingly, an object of the present invention to provide a method for identifying a compound as a GABA<sub>B</sub> receptor agonist or antagonist, said method

5 comprising;

a) administering the compound to a cellular composition of cells expressing a functional GABA<sub>B</sub> receptor, wherein said cells do not normally express the GABA<sub>B</sub> receptor, in the presence of a detectably labeled agonist or antagonist of the GABA<sub>B</sub> receptor; and

10 b) determine the binding of the labeled agonist or antagonist to said cellular composition,

where if the amount of binding of the labeled agonist or antagonist is less in the presence of the test compound, then the compound is a potential agonist respectively antagonist of the GABA<sub>B</sub> receptor.

15 As already specified for the general binding assay above, the binding of the GABA<sub>B</sub> receptor agonist or antagonist is assessed using art-known methods for the study of protein-ligand interactions. The label is generally selected from a radioactive label, a fluorescent label or an enzymatic label, in particular a radiolabel wherein the agonist is selected from the group consisting of <sup>3</sup>H-GABA, <sup>3</sup>H-baclofen and <sup>3</sup>H-3-APPA and the

20 antagonist is selected from the group consisting of <sup>3</sup>H-CGP542626 and <sup>3</sup>H-SCH50911. In a more specific embodiment the aforementioned binding assays are performed on a cellular composition consisting of the membrane fraction of cells according to the invention, in particular on membrane fractions of the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone 20 on August 22, 2003 with the accession

25 number **LMBP XXXX**, using one or more of the aforementioned radiolabeled agonist and/or antagonists.

In a further embodiment the present invention provides a functional assay for

30 identifying compounds that modulate the GABA<sub>B</sub>-receptor activity in the cells according to the invention. Such an assay is conducted using the cells of the present invention, i.e. cotransfected with the human GABA<sub>B</sub>R1a and human GABA<sub>B</sub>R2 subunits. The cells are contacted with at least one reference compound wherein the ability of said compound to modulate the GABA<sub>B</sub>-receptor activity is known. Thereafter, the cells are

35 contacted with a test compound and determined whether said test compound modulates the activity of the GABA<sub>B</sub> receptor compared to the reference compound. A "reference

compound" as used herein refers to a compound that is known to bind and/or to modulate the GABA<sub>B</sub> receptor activity.

5 A compound or a signal that "modulates the activity" of a polypeptide of the invention refers to a compound or a signal that alters the activity of the polypeptide so that it behaves differently in the presence of the compound or signal than in the absence of the compound or signal. Compounds affecting modulation include agonists and antagonists. An agonist of the GABA<sub>B</sub> receptor encompasses a compound such as GABA, baclofen and 3 - APPA which activates GABA<sub>B</sub> receptor function.  
10 Alternatively, an antagonist includes a compound that interferes with GABA<sub>B</sub> receptor function. Typically, the effect of an antagonist is observed as a blocking of agonist-induced receptor activation. Antagonists include competitive as well as non-competitive antagonists. A competitive antagonist (or competitive blocker) interacts with or near the site specific for agonist binding. A non-competitive antagonist or  
15 blocker inactivates the function of the receptor by interacting with a site other than the agonist interaction site.

In one embodiment the present invention provides a method for identifying compounds that have the capability to modulate GABA<sub>B</sub> receptor activity, said method comprising;  
20 a) contacting cells expressing a functional GABA<sub>B</sub> receptor, wherein said cells do not normally express a functional GABA<sub>B</sub> receptor, with at least one reference compound, under conditions permitting the activation of the GABA<sub>B</sub> receptor;  
b) contacting the cells of step a) with a test compound, under conditions permitting the activation of the GABA<sub>B</sub> receptor, and  
25 c) determine whether said test compound modulates the GABA<sub>B</sub> receptor activity compared to the reference compound.

Methods to determine the capability of a compound to modulate the GABA<sub>B</sub> receptor activity are based on the variety of assays available to determine the functional response  
30 of G-protein coupled receptors (see above) and in particular on assays to determine the changes in potassium currents, changes in calcium concentration, changes in cAMP and changes in GTP $\gamma$ S binding. Conditions permitting the activation of the GABA<sub>B</sub> receptor generally known in the art, for example in case of antagonist screening these conditions comprise the presence of a GABA<sub>B</sub> receptor agonist in the assay system.  
35 Typical GABA<sub>B</sub> receptor agonists used in these activity assays are GABA, baclofen or 3-APPA. More particular in the GTP $\gamma$ S assay as outlined herein below, GABA is used

to activate the GABA<sub>B</sub> receptor in order to assess the capability of a test compound to inactivate the GABA<sub>B</sub> receptor protein.

In the aforementioned assay an increase of GTP $\gamma$ S binding in the presence of the test compound is an indication that the compound activates the GABA<sub>B</sub> receptor activity, and accordingly that said test compound is a potential agonist of the GABA<sub>B</sub> receptor protein. A decrease of GTP $\gamma$ S binding in the presence of the test compound is an indication that the compound inactivates the GABA<sub>B</sub> receptor protein and accordingly that said test compound is a potential antagonist of the GABA<sub>B</sub> receptor protein.

Particularly preferred types of assays include binding assays and functional assays which may be performed as follows:

Binding assays

Over-expression of the GABA<sub>B</sub> receptor expressed by the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line of the present invention may be used to produce membrane preparations bearing said receptor (referred to in this section as GABA<sub>B</sub> binding receptor for convenience) for ligand binding studies. These membrane preparations can be used in conventional filter-binding assays (eg. Using Brandel filter assay equipment) or in high throughput Scintillation Proximity type binding assays (SPA and Cytostar-T flashplate technology; Amersham Pharmacia Biotech) to detect binding of radio-labelled GABA<sub>B</sub> ligands (including <sup>3</sup>H-GABA, <sup>3</sup>H-baclofen, <sup>3</sup>H-3-APPA, <sup>3</sup>H-CGP542626, <sup>3</sup>H-SCH50911) and displacement of such radio-ligands by competitors for the binding site. Radioactivity can be measured with Packard Topcount, or similar instrumentation, capable of making rapid measurements from 96-, 384-, 1536- microtitre well formats. SPA/Cytostar-T technology is particularly amenable to high throughput screening and therefore this technology is suitable to use as a screen for compounds able to displace standard ligands.

Another approach to study binding of ligands to GABA<sub>B</sub> binding receptor protein in an environment approximating the native situation makes use of a surface plasmon resonance effect exploited by the Biacore instrument (Biacore). GABA<sub>B</sub> binding receptor in membrane preparations or whole cells could be attached to the biosensor chip of a Biacore and binding of ligands examined in the presence and absence of compounds to identify competitors of the binding site.



Functional assays

Since GABA<sub>B</sub> receptors belong to the family G-protein coupled receptors that are coupled to GIRK (inward rectifying potassium channels), potassium ion flux should result on activation of these receptors. This flux of ions may be measured in real time using a variety of techniques to determine the agonistic or antagonistic effects of particular compounds. Therefore, recombinant GABA<sub>B</sub> binding receptor proteins expressed in the cell lines of the present invention can be characterised using whole cell and single channel electrophysiology to determine the mechanism of action of compounds of interest. Electrophysiological screening, for compounds active at GABA<sub>B</sub> binding receptor proteins, may be performed using conventional electrophysiological techniques and when they become available, novel high throughput methods currently under development.

Given the presynaptic effect of GABA<sub>B</sub> receptor activation on Ca<sup>2+</sup> channels, in an alternative functional screen the modulatory effect of a compound is assessed through the changes in intracellular calcium. Calcium fluxes are measurable using several ion-sensitive fluorescent dyes, including fluo-3, fluo-4, fluo-5N, fura red and other similar probes from suppliers including Molecular Probes. The inhibition of calcium influx as a result of GABA<sub>B</sub> receptor activation can thus be characterised in real time, using fluorometric and fluorescence imaging techniques, including fluorescence microscopy with or without laser confocal methods combined with image analysis algorithms.

Another approach is a high throughput screening assay for compounds active as either agonists or modulators which affect calcium transients. This assay is based around an instrument called a **FLuorescence Imaging Plate Reader (FLIPR®)**, Molecular Devices Corporation). In its most common configuration, it excites and measures fluorescence emitted by fluorescein-based dyes. It uses an argon-ion laser to produce high power excitation at 488 nm of a fluorophore, a system of optics to rapidly scan the over the bottom of a 96-/384-well plate and a sensitive, cooled CCD camera to capture the emitted fluorescence. It also contains a 96-/384-well pipetting head allowing the instrument to deliver solutions of test agents into the wells of a 96-/384-well plate. The FLIPR assay is designed to measure fluorescence signals from populations of cells before, during and after addition of compounds, in real time, from all 96-/384-wells simultaneously. The FLIPR assay may be used to screen for and characterise compounds functionally active at the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line.

A high throughput screening assay, specifically usefull to identify GABA<sub>B</sub> agonists could consist of an arrangement wherein hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cells, are loaded with an appropriate fluorescent dye, incubated with a test compound and  
5 after sufficient time to allow interaction (8 – 24 hours, typically 12-24 hours, in particular 24 hours.) the change in relative fluorescence units measured using an automated fluorescence plate reader such as FLIPR or Ascent Fluoroskan (commercially available from Thermo Labsystems, Brussel, Belgium).

10 In a further embodiment the functional assay is based on the change in GTP $\gamma$ S binding to the GABA<sub>B</sub> binding receptor. In particular using a compotion bindig assay to determine the displacement of radiolabelled GTP $\gamma$ S. In general, this method to identify GABA<sub>B</sub>-receptor agonists comprises preparing a membrane fraction from cells  
15 expressing the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 heterodimer af the present invention, contacting said membrane preparations with the compound to be tested in the presence of radiolabelled GTP $\gamma$ S, under conditions permitting the activation of the GABA<sub>B</sub> receptor, and detecting GTP $\gamma$ S binding to the membrane fraction. An increase in GTP $\gamma$ S binding in the presence of the compound is an indication that the compound  
20 activates the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 receptor. A decrease in GTP $\gamma$ S binding in the presence of the compound is an indication that the compound inactivates the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 receptor. Preferably this GTP $\gamma$ S binding assay is performed on membrane fractions obtained from the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as CHO-  
25 K1 h-GABA-b R1a/R2 clone 20 on August 22, 2003 with the accession number **LMBP XXXX**. Further, the conditions permitting the activation of the GABA<sub>B</sub> receptor comprise the presence of a GABA<sub>B</sub> receptor agonist, such as for example GABA, baclofen and 3-APPA in the assay system. In particular GABA.

This and other functional screening assays will be provided in the examples hereinafter.

30

#### *Method of Treatment*

The present invention also provides the use of a compound identified as a GABA<sub>B</sub> receptor activity modulator, using one of the aforementioned assays, in the manufacture of a medicament for the treatment an indication such as stiff man syndrome,  
35 gastroesophogeal reflux, neuropathic pain, incontinence and treatment of cough and cocaine addiction. In particular for use in the manufacture of a medicament to reduce transient lower esophageal sphincter relaxations (TLESR). It is thus an object of the

present invention to provide a method for the treatment of a warm-blooded animal, for example, a mammal including humans, suffering from an indication such as stiff man syndrome, gastroesophageal reflux, neuropathic pain, incontinence and treatment of cough and cocaine addiction, in particular TLESR.

5        Said method comprising administering to a warm-blooded animal in need thereof an effective amount of a compound identified as a GABA<sub>B</sub> receptor modulator using a method according to the invention. In particular the systemic or topical administration of an effective amount of a compound according to the invention, to warm-blooded animals, including humans.

10        Such agents may be formulated into compositions comprising an agent together with a pharmaceutically acceptable carrier or diluent. The agent may in the form of a physiologically functional derivative, such as an ester or a salt, such as an acid addition salt or basic metal salt, or an N or S oxide. Compositions may be formulated for any  
15        suitable route and means of administration. Pharmaceutically acceptable carriers or diluents include those used in formulations suitable for oral, rectal, nasal, inhalable, topical (including buccal and sublingual), vaginal or parenteral (including subcutaneous, intramuscular, intravenous, intradermal, intrathecal and epidural) administration. The choice of carrier or diluent will of course depend on the proposed  
20        route of administration, which, may depend on the agent and its therapeutic purpose. The formulations may conveniently be presented in unit dosage form and may be prepared by any of the methods well known in the art of pharmacy. Such methods include the step of bringing into association the active ingredient with the carrier which constitutes one or more accessory ingredients. In general the formulations are prepared  
25        by uniformly and intimately bringing into association the active ingredient with liquid carriers or finely divided solid carriers or both, and then, if necessary, shaping the product.

30        For solid compositions, conventional non-toxic solid carriers include, for example, pharmaceutical grades of mannitol, lactose, cellulose, cellulose derivatives, starch, magnesium stearate, sodium saccharin, talcum, glucose, sucrose, magnesium carbonate, and the like may be used. The active compound as defined above may be formulated as suppositories using, for example, polyalkylene glycols, acetylated triglycerides and the like, as the carrier. Liquid pharmaceutically administrable compositions can, for  
35        example, be prepared by dissolving, dispersing, etc, an active compound as defined above and optional pharmaceutical adjuvants in a carrier, such as, for example, water, saline aqueous dextrose, glycerol, ethanol, and the like, to thereby form a solution or

suspension. If desired, the pharmaceutical composition to be administered may also contain minor amounts of non-toxic auxiliary substances such as wetting or emulsifying agents, pH buffering agents and the like, for example, sodium acetate, sorbitan monolaurate, triethanolamine sodium acetate, sorbitan monolaurate, triethanolamine oleate, etc. Actual methods of preparing such dosage forms are known, or will be apparent, to those skilled in this art; for example, see Gennaro et al., Remington's Pharmaceutical Sciences, Mack Publishing Company, Easton, Pennsylvania, 18th Edition, 1990.

- 10 The composition or formulation to be administered will, in any event, contain a quantity of the active compound(s) in an amount effective to alleviate the symptoms of the subject being treated.

15 Dosage forms or compositions containing active ingredient in the range of 0.25 to 95% with the balance made up from non-toxic carrier may be prepared.

20 For oral administration, a pharmaceutically acceptable non-toxic composition is formed by the incorporation of any of the normally employed excipients, such as, for example, pharmaceutical grades of mannitol, lactose, cellulose, cellulose derivatives, sodium crosscarmellose, starch, magnesium stearate, sodium saccharin, talcum, glucose, sucrose, magnesium carbonate, and the like. Such compositions take the form of solutions, suspensions, tablets, pills, capsules, powders, sustained release formulations and the like. Such compositions may contain 1%-95% active ingredient, more preferably 2-50%, most preferably 5-8%.

25 Parenteral administration is generally characterized by injection, either subcutaneously, intramuscularly or intravenously. Injectables can be prepared in conventional forms, either as liquid solutions or suspensions, solid forms suitable for solution or suspension in liquid prior to injection, or as emulsions. Suitable excipients are, for example, water, saline, dextrose, glycerol, ethanol or the like. In addition, if desired, the pharmaceutical compositions to be administered may also contain minor amounts of non-toxic auxiliary substances such as wetting or emulsifying agents, pH buffering agents and the like, such as for example, sodium acetate, sorbitan monolaurate, triethanolamine oleate, triethanolamine sodium acetate, etc.

35 The percentage of active compound contained in such parental compositions is highly dependent on the specific nature thereof, as well as the activity of the compound and

the needs of the subject. However, percentages of active ingredient of 0.1% to 10% in solution are employable, and will be higher if the composition is a solid which will be subsequently diluted to the above percentages. Preferably, the composition will comprise 0.2-2% of the active agent in solution.

5

Throughout this description the terms "standard methods", "standard protocols" and "standard procedures", when used in the context of molecular biology techniques, are to be understood as protocols and procedures found in an ordinary laboratory manual such as: Current Protocols in Molecular Biology, editors F. Ausubel et al., John Wiley and Sons, Inc. 1994, or Sambrook, J., Fritsch, E.F. and Maniatis, T., Molecular Cloning: A laboratory manual, 2nd Ed., Cold Spring Harbor Laboratory Press, Cold Spring Harbor, NY 1989.

10

15 This invention will be better understood by reference to the Experimental Details that follow, but those skilled in the art will readily appreciate that these are only illustrative of the invention as described more fully in the claims that follow thereafter. Additionally, throughout this application, various publications are cited. The disclosure of these publications is hereby incorporated by reference into this application to describe more fully the state of the art to which this invention pertains.

20

## EXPERIMENTAL

### 25 MATERIAL AND METHODS

*Permanent transfection of GABABR1a and GABABR2 in CHO-K1 cells using Lipofectamine PLUS:*

30 CHO-cells were transfected with hGABABR1a/pcDNA3.1. Monoclonal stable R1a-expressing cells were transfected with hGABABR2/pcDNA3.1Hygro+. Selection of clones occurred with 800 µg geneticin + 800 µg hygromycin/ml.

*Membrane preparation:*

35 Butyrate-stimulated (5 mM final) cells were scraped, after a short rinse with PBS, in 50 mM TrisHCl pH7.4 and centrifuged at 23500 g for 10 min. at 4°C. The pellet was homogenised in 5 mM TrisHCl pH 7.4 by Ultra-Turrax (24000 rpm) followed by

centrifugation at 30000 g for 20 min. at 4°C. The resulting pellet was resuspended in 50 mM TrisHCl pH 7.4 and rehomogenised. Protein concentration was determined using the Bradford method.

5 *GTP $\gamma$ 35S activation assay:*

10  $\mu$ g membrane prep was incubated in 250  $\mu$ l in 20 mM Hepes pH 7.4, 100 mM NaCl, 3 mM MgCl<sub>2</sub>, 0.25 nM GTP $\gamma$ 35S, 3  $\mu$ M GDP, 10  $\mu$ g saponin/ml, with or without 1mM GABA (basal activity in absence of baclofen) at 37°C for 20 min. Filtration was carried out onto 96-well GF/B filter plate in Harvester (Packard). Filters were rinsed 6  
10 times with cold 10 mM phosphate buffer pH 7.4, and dried overnight before addition of 30  $\mu$ l Microscint O, and measurement in Topcount (Packard, 1min./well).

*3H-agonist binding:*

30 - 60  $\mu$ g membrane prep was incubated in 50 mM TrisHCl pH 7.4, 2.5 mM CaCl<sub>2</sub>,  
15 10 nM 3H-GABA or 20 nM 3H-baclofen in 500  $\mu$ l at 20°C. Non-specific binding was determined in the presence of 100  $\mu$ M baclofen. After 90 minutes the mixture was transferred onto 96-well GF/B filterplate by Harvester (Packard). Filters were rinsed 6 times with cold 50 mM TrisHCl pH 7.4, 2.5 mM CaCl<sub>2</sub>, and dried overnight before addition of 30  $\mu$ l Microscint O, and measurement in Topcount (Packard, 1min./well).

20

## RESULTS

*GTP $\gamma$ 35S activation assay*

In membranes of stably hGABABR1a-transfected CHO-cells, we measured binding of  
25 the antagonist 3H-CGP54626. hGABABR2 was co-transfected in those R1a-clones with the highest antagonist binding. After subcloning stable clones were obtained showing functional activity in GTP $\gamma$ 35S-binding assay upon stimulation of membranes by GABA, wherein said activity was potentiated in the presence of the positive modulator CGP7930 (Urwyler S., *et al.*, 2001, Molecular Pharmacology 60:963-971)  
30 (fig. 1).

*Agonist Filter Binding Assay*

An agonist filter binding assay has been developed in 96-well GF/B filterplate. The IC<sub>50</sub> of known agonists and antagonists was determined (fig.2). Unexpectedly, in our  
35 hGABABR1a/R2 clone agonist binding was detected with 3H-baclofen as well as with

<sup>3</sup>H-GABA . The K<sub>d</sub> for <sup>3</sup>H-baclofen, <sup>3</sup>H-GABA, and <sup>3</sup>H-CGP54626 was determined in saturation experiments and compared well with published results obtained with tissue preparations (table 1).

Table 1

5

**<sup>3</sup>H-baclofen**

Rat	132 nM	(Hill & Bowery, 1981)
Dog cortex	28 nM	(J&JPRD, 2000)
hGABA <sub>B</sub> R1aR2/CHO	30 nM	(our data, n=2))

10

**<sup>3</sup>H-GABA**

Rat	77 nM	(Hill & Bowery, 1983)
Rat	15-30 nM	(Cross & Horton, 1988)
Pig	26 nM	(Facklam & Bowery, 1993)
Human	20-30 nM	(Cross & Horton, 1988)
hGABA <sub>B</sub> R1aR2/CHO	10-30 nM	(our data, n=6)

15

**<sup>3</sup>H-CGP54626**

Rat	1.5 nM	(Bittiger et al., 1993)
Pig	1.35 nM	(Facklam & Bowery, 1993)
hGABA <sub>B</sub> R1aR2/CHO	1.5 nM	(Green et al., 1993)
hGABA <sub>B</sub> R1aR2/CHO	2.78 nM	(our data, n=1)

20

25 The order of potency for agonists was AMPA > GABA > baclofen, and for antagonists CGP54626 > SCH50911 (fig.2). The obtained IC<sub>50</sub>s were reproducible between different membrane preparations (fig.3)

30 *Agonist centrifugation Binding Assay* In an alternative binding assay the non-bound ligand was separated from the membranes by centrifugation instead of filtration. The assay was performed according to the earlier described filter binding assay, with the difference that the non-bound ligand was separated from the membranes by centrifugation in a microcentrifuge at 12500 rpm for 10 minutes. The supernatant was discarded, the pellet was rinsed with washing buffer and dissolved in 200 µl water. Scintillation fluid was added and the bound <sup>3</sup>H-GABA measured in Topcount (Packhard, 1 min./well).

35

In a saturation assay using increasing concentrations of <sup>3</sup>H-GABA (1 – 400 nM final) I was found that the GABA<sub>B</sub> receptor expressed by the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO

cell line, possess a low and a high affinity agonist binding site. Results of the saturation and scatchard analysis are summarized in Table 2.

Table 2

	Mean (n=5) nM	SD nM
<b>Bmax 1</b>	0.19	0.05
<b>Kd 1</b>	9.4	3.1
<b>Bmax 2</b>	0.76	0.24
<b>Kd 2</b>	401	224

## DISCUSSION

To our knowledge, no earlier reports were made in literature of recombinant hGABA<sub>B</sub> receptor, showing agonist binding with a high and low affinity binding site in a filter binding assay. An HTS agonist filter binding screen has been developed using 3H-GABA. We found reproducible K<sub>i</sub> values for known agonists and antagonists, independent of the membrane preparation.

It has in addition been demonstrated that the recombinant GABA<sub>B</sub> receptor has two agonist binding sites. One high affinity and one low affinity binding site. It is to be expected that high affinity agonists of the GABA<sub>B</sub> receptor will elicit a different response compared to the low affinity agonists. Hence, the cell line of the present invention not only allows to identify GABA<sub>B</sub> receptor agonists, but also provides a useful tool to characterize the nature of the compound receptor interaction.



SEQUENCE LISTING

5 <110> Janssen Pharmaceutica N.V.

<120> CHIMERIC GABA<sub>A</sub> RECEPTOR

10 <130> PRD 2108p

<160> 4

15 <210> 1  
<211> 2886  
<212> DNA  
<213> Homo sapiens

20 <220>  
<221> CDS  
<222> (1)..(2886)  
<223>

25 <400> 1  
atgttgctgc tgctgtact ggcggcactc ttcctccgcc ccccgggcgc gggcggggcg 60  
cagacccccca acgccacctc agaagggtgc cagatcatat acccgccctg ggaagggggc 120  
atcagggtacc ggggcctgac tcgggaccag gtgaaggcta tcaacttcct gccagtggac 180  
30 tatgagattg agtatgtgtg ccgggggggag cgcgagggtg tggggcccaa ggtccgcaag 240  
tgcttgccca acggctcctg gacagatatg gacacacca gccgctgtgt ccgaatctgc 300  
tccaagtctt atttgaccct ggaaaatggg aagggtttcc tgacgggtgg ggacctccca 360  
gctctggacg gagcccggtt ggatttccgg tgtgacctcg acttccatct ggtgggcagc 420  
tcccggagca tctgtagtca gggccagtgg agcaccacca agccccactg ccagggtgaat 480  
35 cgaacgccac actcagaacg gcgcgcagtg tacatcgggg cactgtttcc catgagcggg 540  
ggctggccag ggggcccaggc ctgccagccc gcggtggaga tggcgctgga ggacgtgaat 600  
agccgcaggg acatcctgcc ggactatgag ctcaagctca tccaccacga cagcaagtgt 660  
gatccaggcc aagccaccaa gtacctatat gagctgtctt acaacgacct tatcaagatc 720  
40 atccttatgc ctggctgcag ctctgtctcc acgctgggtg ctgaggctgc taggatgtgg 780  
aacctcattg tgctttccta tggatccagc tcaccagccc tgtcaaaccg gcagcgtttc 840  
cccactttct tccgaacgca cccatccagc acactccaca accctaccog cgtgaaactc 900  
tttgaagaag ggggctggaa gaagattgct accatccagc agaccactga ggtcttctact 960  
tcgactctgg acgacctgga ggaacgagtg aaggaggctg gaattgagat tactttccgc 1020  
45 cagagtttct tctcagatcc agctgtgccc gtcaaaaacc tgaagcgcca ggatgccga 1080  
atcatcgtgg gacttttcta tgagactgaa gcccggaag ttttttgtga ggtgtacaag 1140  
gagcgtctct ttgggaagaa gtacgtctgg ttcctcattg ggtgggtatg tgacaattgg 1200  
ttcaagatct acgaccttcc tatcaactgc acagtggatg agatgactga ggcgggtggg 1260  
ggccacatca caactgagat tgtcatgctg aatcctgcca ataccgcag catttccaac 1320  
50 atgacatccc aggaatttgt ggagaaacta accaagcgac tgaaaagaca ccctgaggag 1380  
acaggaggct tccaggaggc accgctggcc tatgatgcca tctgggcctt ggcactggcc 1440  
ctgaacaaga catctggagg aggcggccgt tctgggtgtg gcctggagga cttcaactac 1500  
aacaaccaga ccattaccga ccaaatctac cgggcaatga actcttcgtc ctttgagggt 1560  
gtctctggcc atgtgggtgt tgatgccagc ggctctcgga tggcatggac gcttatcgag 1620  
55 cagcttcagg gtggcagcta caagaagatt ggctactatg acagcaccaa ggatgatctt 1680  
tcctgggtcca aaacagataa atggattgga ggggtccccc cagctgacca gaccctgggtc 1740  
atcaagacat tccgcttccct gtcacagaaa ctctttatct ccgtctcagt tctctccagc 1800  
ctggggcattg tctcagctgt tgtctgtctg tcccttaaca tctacaactc acatgtccgt 1860  
tatatccaga actcagacg caacctgaac aacctgactg ctgtgggctg ctcactgggt 1920  
60 ttagctgtctg tcttccccct ggggctcgat gggtaccaca ttgggaggaa ccagtttcct 1980  
ttcgtctgcc agggcccgct ctggctcctg ggcctgggct ttagtctggg ctacgggtcc 2040  
atgttcacca agatattgggt ggtccacacg gtcttcacaa agaaggaaga aaagaaggag 2100  
tgaggagaaga ctctggaacc ctggaagctg tatgccacag tgggcctgct ggtgggcatg 2160  
gatgtcctca ctctcgccat ctggcagatc gtggaccctc tgcaccggac cattgagaca 2220  
75 tttgccaagg aggaacctaa ggaagatatt gacgtctcta ttctgcccc gctggagcat 2280  
tgcagctcca ggaagatgaa tacatggctt ggcattttct atggttaciaa ggggctgctg 2340  
ctgctgtctgg gaatcttccct tgccttatgag accaagagtg tgtccactga gaagatcaat 2400  
gatcaccggg ctgtgggcat ggcttatctac aatgtggcag tcctgtgcct catcactgct 2460  
cctgtcacca tgattctgtc cagccagcag gatgcagcct ttgcctttgc ctctcttgc 2520  
atagttttct cctcctatat cactcttgtt gtgctctttg tgcccaagat gcgcaggctg 2580

5	atcaccccgag	gggaatggca	gtcggaggcg	caggacacca	tgaagacagg	gtcatcgacc	2640
	aacaacaacg	aggaggagaa	gtcccggctg	ttggagaagg	agaaccgtga	actggaaaag	2700
	atcattgctg	agaaagagga	gcgtgtctct	gaactgcgcc	atcaactcca	gtctcggcag	2760
	cagctccgct	cccgccgcca	ccccaccgaca	cccccagaac	cctctggggg	cctgcccagg	2820
	ggaccacctg	agccccccga	ccggcttagc	tgtgatggga	gtcagagtga	tttgctttat	2880
	aagtga						2886
10	<210>	2					
	<211>	961					
	<212>	PRT					
	<213>	Homo sapiens					
15	<400>	2					
	MLLLLLLAPL	FLRPPGAGGA	QTPNATSEGC	QIIHPPWEGG	IRYRGLTRDQ	VKAINFLPVD	60
	YEIEYVCRGE	REVVGPKVRK	CLANGSWTDM	DTPSRCVRIK	SKSYLTLENG	KVFLTGGDLP	120
	ALDGARVDFR	CDFDFHLVGS	SRSICSQGW	STPKPHCQVN	RTPHSERRAV	YIGALFPMSG	180
20	GWPGGQACQP	AVEMALEDVN	SRRDILPDYE	LKLIHDSKC	DPGQATKYLY	ELLYNDPIKI	240
	ILMPGCSSVS	TLVAEAARMW	NLIVLSYGSS	SPALSNRQRF	PTFFRTHPSA	TLHNPTRVKL	300
	FEKWGWKKIA	TIQQTTEVFT	STLDDLEERV	KEAGIEITFR	QSFFSDPAVP	VKNLKRQDAR	360
	IIVGLFYETE	ARKVFCEVYK	ERLFGKKYVW	FLIGWYADNW	FKIYDPSINC	TVDEMTEAVE	420
25	GHIITTEIVML	NPANTRISIS	MTSQEFVEKL	TKRLKRHPPE	TGGFQEAFLA	YDAIWALALA	480
	LNKTSGGGGR	SGVRLEDFNY	NNQTITDQIY	RAMNSSSFEG	VSGHVVDAS	GSRMAWTIE	540
	QLQGSYKIK	GYDSTKDDL	SWSKTDKWIG	GSPPADQTLV	IKTFRFLSQK	LFISVSVLSS	600
	LGIVLAVVCL	SFNIYNSHVR	YIQNSQPNLN	NLTAVGCSLA	LAAVFPLGLD	GYHIGRNQFP	660
30	FVCQARLWLL	GLGFSLGYGS	MFTKIWWVHT	VPTKKEEKKE	WRKTLEPWKL	YATVGLLVGM	720
	DVLTALAIWQI	VDPLHRTIET	FAKEEPKEDI	DVSILPQLEH	CSSRKMNTWL	GIFYGYKGLL	780
	LLLGIFLAYE	TKSVSTEKIN	DHRAVGMAIY	NVAVLCILITA	PVTMILSSQK	DAAFASFALA	840
	IVFSSYITLV	VLFPVKMRL	ITRGEWQSEA	QDTMTKGSST	NNNEEEKSRL	LEKENRELEK	900
	IIAEKEERVS	ELRHQLQSRQ	QLRSRRHPPT	PPEFSGGLPR	GPPEPPDRLS	CDGSRVHLLY	960
	K						961
35	<210>	3					
	<211>	2823					
	<212>	DNA					
	<213>	Homo sapiens					
40	<220>						
	<221>	CDS					
	<222>	(1) .. (2823)					
	<223>						
45	<400>	3					
	atggcttccc	cgcgagactc	cgggcagccc	ggggcgccgc	cgccgcgcgc	accgcccgcc	60
	gcgcgcctgc	tactgtact	gctgctgccg	ctgctgctgc	ctctggcgcc	cggggcctgg	120
	ggctggcgcc	ggggcgcccc	ccggcgcccg	cccagcagcc	cgccgctctc	catcatgggc	180
50	ctcatgccgc	tcaccaagga	ggtggccaag	ggcagcatcg	ggcgcggtgt	gtccccgc	240
	gtggaaactgg	ccatcgagca	gatccgcaac	gagtcactcc	tgccccccta	cttctcgac	300
	ctgcggtctc	atgacacgga	gtgcgacaac	gcaaaagggt	tgaaagcctt	ctacgatgca	360
	ataaaaatcag	ggccgaacca	cttgatgggtg	tttggaggcg	tctgtccatc	cgtcacatcc	420
55	atcattgcag	agtcctccca	aggctggaat	ctgggtgcagc	tttcttttgc	tgcaaccacg	480
	cctgttctag	ccgataagaa	aaaataccct	tatttcttct	ggaccgtccc	atcagacaat	540
	gcggtgaatc	cagccattct	gaagtgtgctc	aagcactacc	agtgggaagcg	cgtgggcacg	600
	ctgacgcaag	acgttcagag	gttctctgag	gtgcggaatg	acctgactgg	agttctgtat	660
60	ggcgaggaca	ttgagatttc	agacaccgag	agcttctcca	acgatccctg	taccagtgtc	720
	aaaaagctga	aggggaatga	tgtgcggatc	atccttggcc	agtttgacca	gaatatggca	780
	gcaaaagtgt	tctgttgtgc	atacaggagg	aacatgtatg	gtagtaaata	tcagtggatc	840
	attccgggct	ggtacgagcc	ttcttgggtg	gagcaggtgc	acacggaagc	caactcatcc	900
65	cgctgcctcc	ggaagaatct	gcttgcctgc	atggagggct	acattggcgt	ggatttcgag	960
	cccctgagct	ccaagcagat	caagaccatc	tcaggaaaga	ctccacagca	gtatgagaga	1020
	gagtacaaca	acaagcggtc	aggcgtgggg	ccagcaagt	tccacgggta	cgcctacgat	1080
	ggcatctggg	tcactgcgca	gacactgcag	agggccatgg	agacactgca	tgccagcagc	1140
70	cgccaccagc	ggatccagga	cttcaactac	acggaccaca	cgctgggcag	gatcatcctc	1200
	aatgccatga	acgagaccaa	cttctctcgg	gtcaccgggtc	aagttgtatt	ccggaatggg	1260
	gagagaatgg	ggaccattaa	atttactcaa	tttcaagaca	gcaggagggt	gaaggtggga	1320
	gagtacaacg	ctgtggccga	cacactggag	atcatcaatg	acaccatcag	gttccaagga	1380
75	tccgaaccac	caaaagacaa	gaccatcatc	ctggagcagc	tgccgaagat	ctccctacct	1440
	ctctacagca	tcctctctgc	cctcaccatc	ctcgggatga	tcattggccag	tgcttttctc	1500
	ttcttcaaca	tcaagaaccg	gaatcagaag	ctcataaaga	tgctcagatcc	atacatgaac	1560
	aaccttatca	tccttggagg	gatgctctcc	tatgcttcca	tatttctctt	tgcccttgat	1620
	ggatcctttg	tctctgaaaa	gacctttgaa	acactttgca	ccgtcaggac	ctggattctc	1680

	accgtgggct	acacgaccgc	ttttggggcc	atgtttgcaa	agacctggag	agtccacgcc	1740
	atcttcaaaa	atgtgaaaat	gaagaagaag	atcatcaagg	accagaaact	gcttgtgatc	1800
	gtggggggca	tgctgctgat	cgacctgtgt	atcctgatct	gctggcaggc	tgtggacccc	1860
5	ctgccaagga	cagtgagaga	gtacagcatg	gagccggacc	cagcaggacg	ggatatctcc	1920
	atccgccctc	tcctggagca	ctgtgagaac	acccatatga	ccatctggct	tggcatcgct	1980
	tatgcctaca	agggacttct	catgttggtc	ggttggttct	tagcttggga	gacccgcaac	2040
	gtcagcatcc	ccgactcaa	cgacagcaag	tacatcgga	tgagtgtcta	caacgtgggg	2100
	atcatgtgca	tcacggggc	cgctgtctcc	ttcctgaccc	gggaccagcc	caatgtgcag	2160
10	ttctgcatcg	tggtctgtgt	catcatcttc	tgcagacca	tcacctctg	cctgggtattc	2220
	gtgcccgaagc	tcacaccct	gagaacaaac	ccagatgcag	caacgcagaa	caggcgattc	2280
	cagttcactc	agaatcagaa	gaaagaagat	tctaaaacgt	ccacctcgt	caccagtgtg	2340
	aaccaagcca	gcacatccc	cctggagggc	ctacagtcag	aaaaccatcg	cctgcgaatg	2400
	aagatcacag	agctggataa	agacttggaa	gaggtcacca	tgacagctga	ggacacacca	2460
	gaaaagacca	cctacattaa	acagaaccac	taccaagagc	tcaatgacat	cctcaacctg	2520
15	ggaaacttca	ctgagagcac	agatggagga	aaggccattt	taaaaaatca	cctcgatcaa	2580
	aatccccagc	tacagtggaa	cacaacagag	ccctctcgaa	catgcaaaga	tcctatagaa	2640
	gatataaact	ctccagaaca	catccagcgt	cggctgtccc	tccagctccc	catcctccac	2700
	cacgcctacc	tcccatccat	cggaggcgtg	gacgccagct	gtgtcagccc	ctgcgtcagc	2760
20	cccaccgcca	gccccgcca	cagacatgtg	ccacctcct	tccgagtcac	ggtctcgggc	2820
	ctg						2823

	<210>	4
	<211>	941
25	<212>	PRT
	<213>	Homo sapiens

	<400>	4						
30	MASPRSSGQP	GPPPPPPPPP	ARLLLLLLLL	LLLPLAPGAW	GWARGAPRPP	PSSPPLSIMG		60
	LMPLTKEVAK	GSIGRGVLP	VELAIEQIRN	ESLLRPYFLD	LRLYDTECDN	AKGLKAFYDA		120
	IKYGNHLMV	FGGVCPSVTS	IIAESLQGN	LVQLSFAATT	PVLADKKKYP	YFFRTVPSDN		180
	AVNPAILKLL	KHYQWKRVGT	LTQDVQRFSE	VRNDLTGVLY	GEDIEISDTE	SFSNDPCTSV		240
	KKLKGNVRI	ILQFDQNM	AKVFCCAYEE	NMYGSKYQWI	IPGWYEPSWW	EQVHTEANSS		300
35	RCLRKNLLAA	MEGYIGVDFE	PLSSQIKITI	SGKTPQQYER	EYNNKRSVG	PSKFHGYAYD		360
	GIWVIAKTLO	RAMETLHASS	RHQRIQDFNY	TDHTLGRIIL	NAMNETNFFG	VTGQVVRNG		420
	ERMGTIKFTQ	FQDSREVKVG	EYNAVADTLE	IINDTIRFQG	SEPPKDKTII	LEQLRKISLP		480
	LYSILSALTI	LGMIMASAF	FFNIKRNQK	LIKMSPPYMN	NLIILGGMLS	YASIFLFLD		540
	GSFVSEKTFE	TLCTVRTWIL	TVGYTTAFGA	MFAKTWRVHA	IFKNVKMKKK	IIKDQKLLVI		600
	VGGMLLIDLC	ILICWQAVDP	LRRTVEKYSM	EPDPAGRDIS	IRPLEHCEN	THMTIWLGI		660
40	YAYKGLMLF	GCFLAWETRN	VSIPALNDSK	YIGMSVYNVG	IMCIIGAAS	FLTRDQPNVQ		720
	FCIVALVIIF	CSTITLCLVF	VPKLITLRTN	PDAATQNRFF	QFTQNQKKED	SKTSTSVTSV		780
	NQASTSRLEG	LQSENHRLRM	KITELDKDLE	EVTMLQDTP	EKTTYIKQNH	YQELNDILNL		840
	GNFTSTDDG	KAILKNHLDQ	NPQLQWNTTE	PSRTCKDPIE	DINSPEHIQR	RLSLQLPILH		900
45	HAYLPSIGGV	DASCVSPCVS	PTASPRHRHV	PPSFRVMVSG	L			941

WHAT IS CLAIMED IS:

1. An isolated GABA<sub>B</sub> receptor protein comprising at least one GABA<sub>B</sub>R1a subunit and at least one GABA<sub>B</sub>R2 subunit, characterized in that said GABA<sub>B</sub> receptor has one high affinity agonist binding site and one low affinity agonist binding site.  
5
2. The GABA<sub>B</sub> receptor protein according to claim 1 wherein the GABA<sub>B</sub>R1a subunit is encoded by the oligonucleotide sequence consisting of SEQ ID No.1 and the GABA<sub>B</sub>R2 subunit is encoded by the oligonucleotide sequence consisting of SEQ ID N0.3.  
10
3. The GABA<sub>B</sub> receptor protein according to claims 1 or 2 wherein said receptor protein is expressed by the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone 20 on August 22, 2003 with the accession number LMBP  
15
4. Use of the GABA<sub>B</sub> receptor protein according to any one of claims 1 to 3 in a method to identify GABA<sub>B</sub> receptor agonists or antagonists.  
20
5. The hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone on August 22, 2003 with the accession number LMBP .
- 25 6. A method to identify whether a test compound binds to a GABA<sub>B</sub> receptor protein according to any one of claims 1 to 3, and is thus a potential agonist or antagonist of the GABA<sub>B</sub> receptor, said method comprising:  
a) contacting cells expressing a functional GABA<sub>B</sub> receptor, wherein such cells do not normally express the GABA<sub>B</sub> receptor, with the test compound in the presence and absence of a compound known to bind to the GABA<sub>B</sub> receptor, and  
30

b) determine the binding of the test compound to the GABA<sub>B</sub> receptor using the compound known to bind to the GABA<sub>B</sub> receptor as a reference.

5 7. A method according to claim 6, wherein the compound known to bind to the GABA<sub>B</sub> receptor is detectably labeled, and wherein said label is used to determine the binding of the test compound to the GABA<sub>B</sub> receptor.

10 8. A method according to claim 7 wherein the compound known to bind to the GABA<sub>B</sub> receptor is selected from the group consisting of <sup>3</sup>H-GABA, <sup>3</sup>H-baclofen, <sup>3</sup>H-3-APPA, <sup>3</sup>H-CGP542626 and <sup>3</sup>H-SCH50911.

15 9. A method to identify GABA<sub>B</sub> receptor agonists said method comprising,  
a) exposing cells expressing a functional GABA<sub>B</sub> receptor, wherein such cells do not normally express the GABA<sub>B</sub> receptor, to a labeled agonist of GABA<sub>B</sub> in the presence and absence of the test compound, and  
b) determine the binding of the labeled agonist to said cells,  
where if the amount of binding of the labeled agonist is less in the presence of the test compound, then the compound is a potential agonist of the GABA<sub>B</sub> receptor.

20 10. A method according to claim 10 wherein the labeled agonist is selected from the group consisting of <sup>3</sup>H-GABA, <sup>3</sup>H-baclofen and <sup>3</sup>H-3-APPA.

25 11. A method to identify GABA<sub>B</sub> receptor antagonists said method comprising,  
a) exposing cells expressing a functional GABA<sub>B</sub> receptor, wherein such cells do not normally express the GABA<sub>B</sub> receptor, to a labeled antagonist of GABA<sub>B</sub> in the presence and absence of the test compound, and  
b) determine the binding of the labeled antagonist to said cells,

where if the amount of binding of the labeled antagonist is less in the presence of the test compound, then the compound is a potential antagonist of the GABA<sub>B</sub> receptor.

- 5    12.    A method according to claim 10 wherein the labeled antagonist is selected from the group consisting of <sup>3</sup>H-CGP542626 and <sup>3</sup>H-SCH50911.
13.    A method for identifying a compound as a GABA<sub>B</sub> receptor agonist, said method comprising;
- 10    a) administering the compound to a cellular composition of the cells according to claim 5, in the presence of a detectably labeled GABA<sub>B</sub> receptor agonist; and
- 15    b) determine the binding of the labeled agonist to said cellular composition, where if the amount of binding of the labeled agonist is less in the presence of the test compound, then the compound is a potential agonist of the GABA<sub>B</sub> receptor.
14.    A method according to claim 13 wherein the cellular composition consists of a membrane fraction of the cells according to claim 5.
- 20    15.    A method according to claims 13 or 14 wherein the labelled agonist is selected from the group consisting of <sup>3</sup>H-GABA, <sup>3</sup>H-baclofen and <sup>3</sup>H-3-APPA.
- 25    16.    A method for identifying a compound as a GABA<sub>B</sub> receptor antagonist, said method comprising;
- 25    a) administering the compound to a cellular composition of the cells according to claim 5, in the presence of a detectably labeled GABA<sub>B</sub> receptor antagonist; and
- 25    b) determine the binding of the labeled antagonist to said cellular composition,

where if the amount of binding of the labeled antagonist is less in the presence of the test compound, then the compound is a potential antagonist of the GABA<sub>B</sub> receptor.

- 5 17. A method according to claim 16 wherein the cellular composition consists of a membrane fraction of the cells according to claim 5.
18. A method according to claims 16 or 17 wherein the labeled antagonist is selected from the group consisting of <sup>3</sup>H-CGP542626 and <sup>3</sup>H-SCH50911.
- 10 19. A method for identifying compounds that have the capability to modulate GABA<sub>B</sub> receptor activity, said method comprising;
  - 15 a) contacting cells expressing a functional GABA<sub>B</sub> receptor, wherein said cells do not normally express a functional GABA<sub>B</sub> receptor, with at least one reference compound, under conditions permitting the activation of the GABA<sub>B</sub> receptor;
  - b) contacting the cells of step a) with a test compound, under conditions permitting the activation of the GABA<sub>B</sub> receptor, and
  - 20 c) determine whether said test compound modulates the GABA<sub>B</sub> receptor activity compared to the reference compound.
20. A method according to claim 19 wherein the capability of the test compound to modulate the GABA<sub>B</sub> receptor activity is determined using one or more of the functional responses selected from the group consisting of changes in potassium  
25 currents, changes in calcium concentration, changes in cAMP and changes in GTPγS binding
21. A method for identifying compounds that have the capability to modulate GABA<sub>B</sub> receptor activity, said method comprising;
  - 30 a) contacting a membrane fraction of the cells according to claim 5, with the compound to be tested in the presence of radiolabeled GTPγS, under conditions permitting the activation of the GABA<sub>B</sub> receptor; and
  - b) determine GTPγS binding to the membrane fraction,  
where an increase in GTPγS binding in the presence of the compound is an  
35 indicaton that the compound activates the GABA<sub>B</sub> receptor activity.

22. A method for identifying compounds that have the capability to modulate GABA<sub>B</sub> receptor activity, said method comprising;
- 5 a) contacting a membrane fraction of the cells according to claim 5, with the compound to be tested in the presence of radiolabelled GTP $\gamma$ S, under conditions permitting the activation of the GABA<sub>B</sub> receptor; and
- b) determine GTP $\gamma$ S binding to the membrane fraction, where an decrease in GTP $\gamma$ S binding in the presence of the compound is an indicaton that the compound inactivates the GABA<sub>B</sub> receptor activity.
- 10 23. A method according to claims 21 or 22 wherein the conditions permitting the activation of the GABA<sub>B</sub> receptor comprise the presence of a GABA<sub>B</sub> receptor agonist.
- 15 24. A method according to claim 23 wherein the GABA<sub>B</sub> receptor agonist is selected from the group consisting of GABA, baclofen and 3-APPA.



## ABSTRACT

5

CHIMERIC GABA<sub>B</sub> RECEPTOR

10 The present invention provides an isolated GABA<sub>B</sub> receptor protein comprising at least one GABA<sub>B</sub>R1a subunit and at least one GABA<sub>B</sub>R2a subunit, characterized in that said GABA<sub>B</sub> receptor has one high affinity agonist binding site and one low affinity agonist binding site. In particular the isolated recombinant GABA<sub>B</sub> receptor protein expressed by the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as  
15 CHO-K1 h-GABA-b R1a/R2 clone on August 22, 2003 with the accession number **LMBP XXXX**. It is thus an object of the present invention to provide the hGABA<sub>B</sub>R1a/GABA<sub>B</sub>R2 CHO cell line deposited at the Belgian Coordinated Collections of Microorganisms (BCCM) as CHO-K1 h-GABA-b R1a/R2 clone on August 22, 2003 with the accession number **LMBP XXXX**.

20 The invention also provides the use of the aforementioned cell line in a method to identify GABA<sub>B</sub> receptor agonists using a functional or a binding assay. In particular in a radioligand-binding assay comprising the use of radiolabeled agonists such as for example <sup>3</sup>H-GABA or <sup>3</sup>H-baclofen.

25 In a particular embodiment the present invention provides the use of the aforementioned GABA<sub>B</sub> receptor in a method to identify a high affinity GABA<sub>B</sub> receptor agonist using a functional or a binding assay. In particular in a radioligand-binding assay comprising the use of radiolabeled agonists such as for example <sup>3</sup>H-GABA or <sup>3</sup>H-baclofen. Alternatively, the aforementioned binding assays are performed on cellular extracts, in particular cellular membrane  
30 preparations of the aforementioned cells.

Fig 1

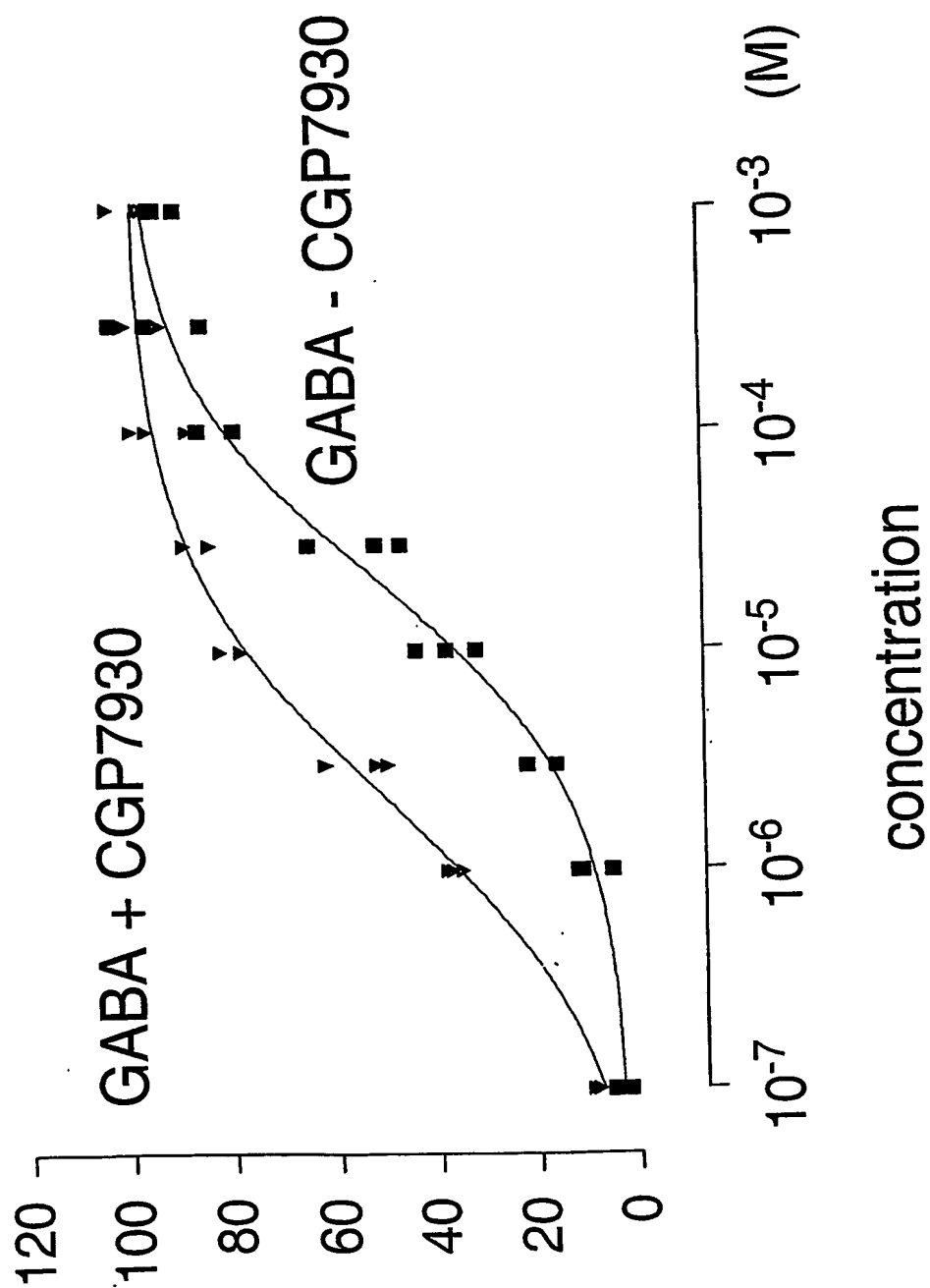


Fig 2

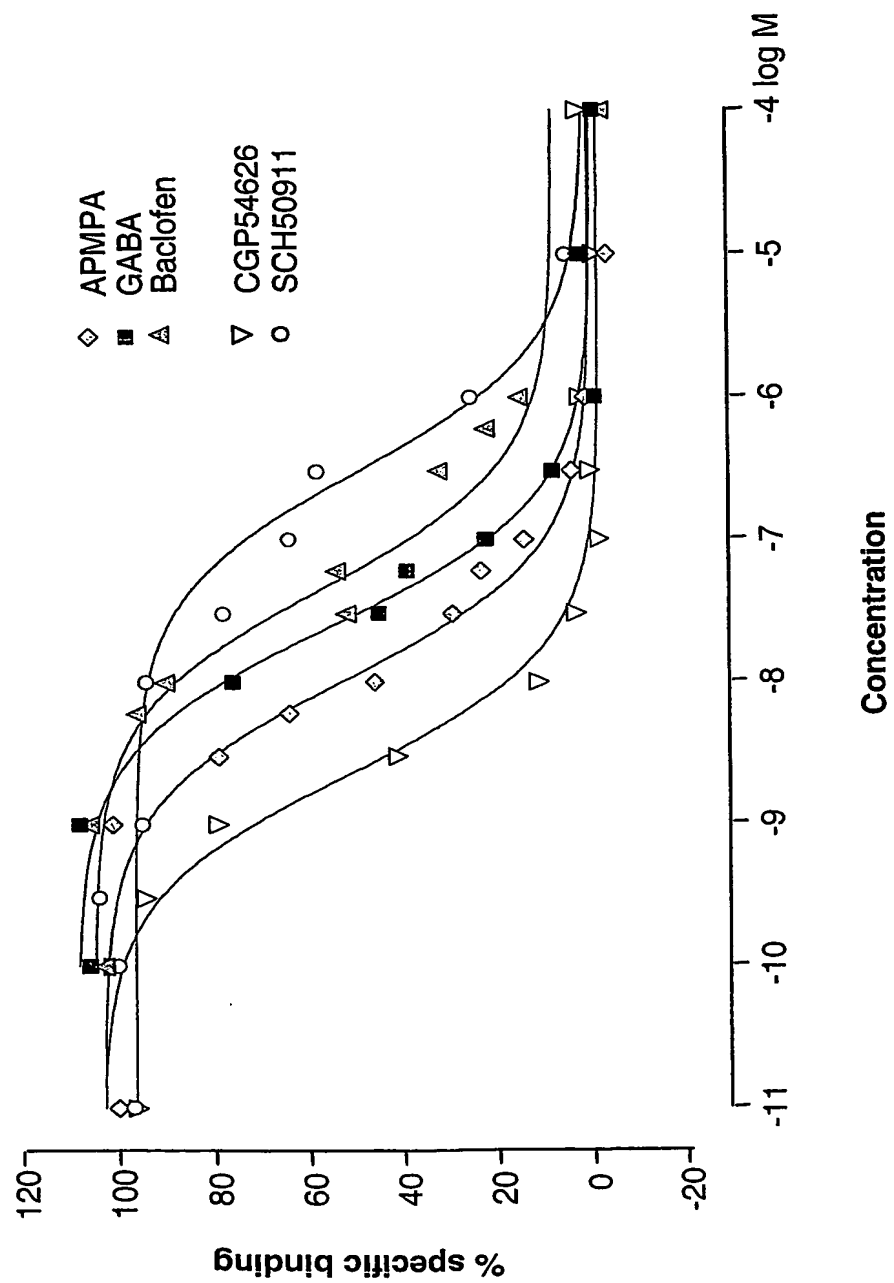


Fig 3

